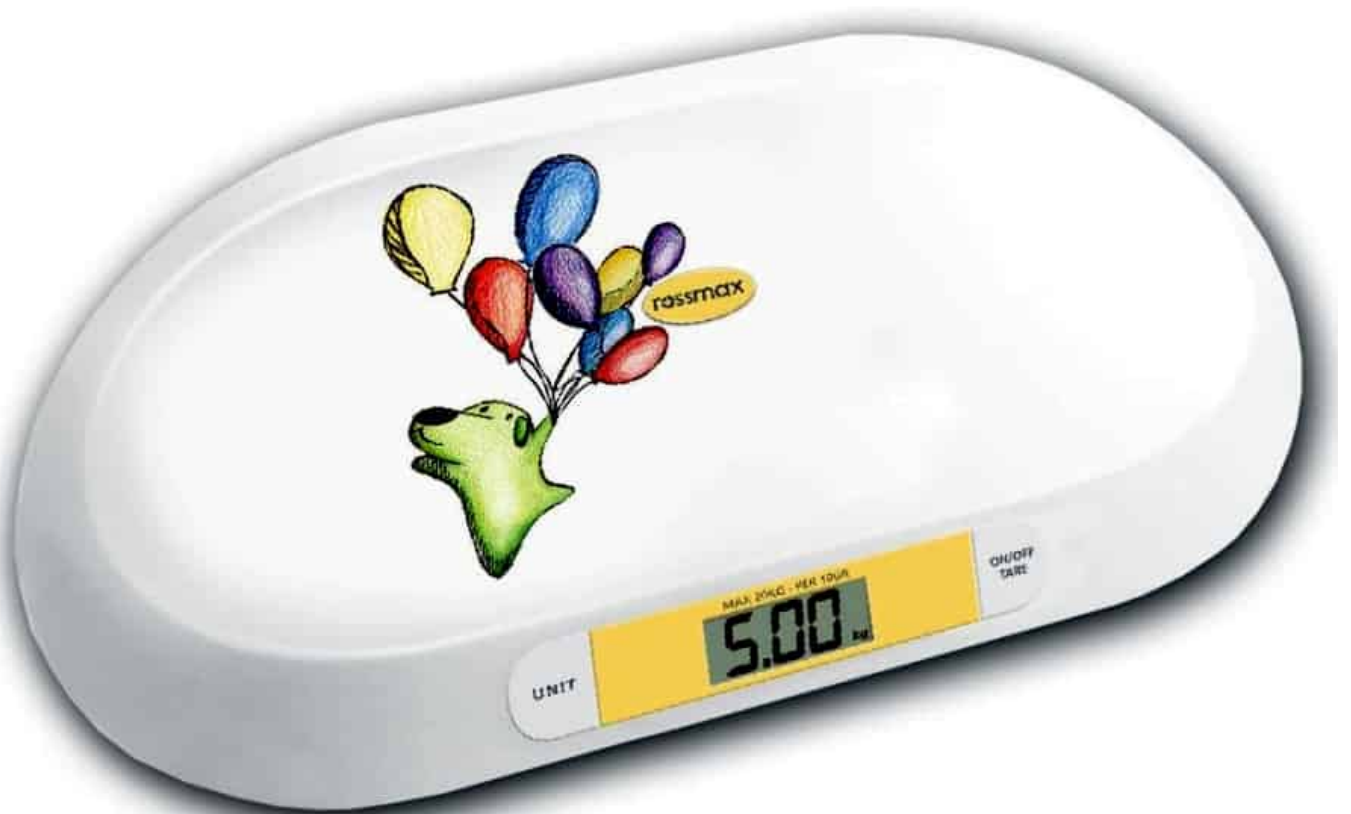




The Rossmax WE300 Qutie Digital Baby Scale (VM-SCAWE300) is a precision infant weighing scale with a 20kg capacity at 5g graduation, designed for accurate weight monitoring of newborns, infants, and young children in clinical and home settings. It features a discretionary hold (HALT) function for active babies, tare function, switchable kg/lb/oz units, auto shut-off, and a low battery indicator, on a

flat 320 × 520mm platform.

Precision Weighing for Newborns and Infants

The Rossmax WE300 Qutie is a Swiss-designed digital baby scale purpose-built for accurate weight monitoring of newborns, infants, and toddlers in paediatric clinics, neonatal units, maternity wards, child and maternal health services, and home settings. With a maximum capacity of 20kg and a fine 5g graduation, the WE300 delivers the measurement precision that infant weight monitoring demands — where changes as small as a few grams carry meaningful clinical significance for feeding assessment, hydration status, and growth velocity tracking. The generously sized 320 × 520mm flat platform provides a safe, stable lying surface from newborn through to toddler age, while the low 160mm overall height keeps the platform close to the support surface for safe placement and removal of the infant.

Practical Clinical Features

The WE300's discretionary hold (HALT) function is its most clinically significant feature. When activated, the scale analyses weight readings of a moving infant and automatically locks the display once a stable measurement is determined, eliminating the challenge of reading an accurate weight from an active or unsettled baby. The tare function allows the weight of a nappy, towel, or light clothing placed on the platform before weighing to be zeroed out, ensuring only the infant's net weight is recorded — a practical necessity in cold environments where full undressing is not appropriate. Measurement units switch between kilograms, pounds, and ounces at the touch of a button, accommodating different clinical documentation standards and caregiver preferences. An overload indicator activates when weight placed on the platform exceeds the 20kg maximum.

Power Management and Reliability

The WE300 is powered by a single CR2032 battery, with an automatic shut-off function conserving battery life between uses. The low battery indicator provides advance warning before battery failure, preventing unexpected interruption during a weighing session — an important consideration where

the scale is used repeatedly throughout a clinical day. The WE300 carries Rossmax's Swiss design accreditation and is backed by a two-year manufacturer's warranty. Available in South Africa from Vital Medical.

TECHNICAL SPECIFICATIONS TABLE

Specification	Detail
Brand	Rossmax
Model	WE300 Qutie
SKU	VM-SCAWE300
Type	Digital baby scale
Maximum Capacity	20kg
Minimum Capacity	100g
Graduation	5g
Units	kg / lb / oz (switchable)
Hold Function	Discretional (HALT)
Tare Function	Yes
Overload Indicator	Yes
Auto Shut-Off	Yes
Low Battery Indicator	Yes
Display	LCD
Platform Dimensions	320(L) × 520(W)mm
Overall Dimensions	320(L) × 525(W) × 160(H)mm
Battery	1 × CR2032
Warranty	2 years

SKU: VM-SCAWE300